



# Ultrasonic Baths

Version 1.0

# Ultrasonic Baths

Our range of ultrasonic cleaning baths is designed for laboratory and industrial use, using high-frequency cavitation and optional heating to efficiently remove contaminants from instruments, components, and delicate items.

## FEATURES + BENEFITS

- Programmed with degassing modes to maintain performance of cavitation
- Equipped with industrial-level transducers to ensure durability
- Lower frequency, at 35kHz, performs extraordinarily well at cleaning small and hard-to-reach crevices
- Tank and housing are made from stainless steel which makes them durable and easy to clean
- Heating units available



ROCKER

## OPERATION



### Applications:

- Solvent degassing: HPLC mobile phase
- Vessel and tool cleaning: cuvette, tubes, glassware, optical lenses and medical tools
- Component cleaning: PCB, electronics and hardware parts

# Ultrasonic Baths

## TECHNICAL SPECIFICATIONS

SPECIFICATIONS	3L ULTRASONIC BATH		6L ULTRASONIC BATH		10L ULTRASONIC BATH		20L ULTRASONIC BATH	
PRODUCT CODE	188203-22	188204-22 (with heating)	188206-22	188207-22 (with heating)	188210-22	188211-22 (with heating)	188220-22	188221-22 (with heating)
TANK CAPACITY	3L		6L		10L		20L	
TANK AND HOUSING MATERIAL	SS304 Stainless Steel							
FREQUENCY	35kHz							
TIMER	1 - 99 minutes							
HEATING	-	R.T. to 60°C	-	R.T. to 60°C	-	R.T. to 60°C	-	R.T. to 60°C
DEGAS	Yes							
DRAIN	-				Yes			
WEIGHT	4kg	4.1kg	5.5kg	5.7kg	7.5kg	7.7kg	14.5kg	15kg
TANK DIMENSION (LxWxH)	24x14x10cm		30x15x15cm		30x24x15cm		50x30x15cm	
OVERALL SIZE (LxWxH)	27x16x24cm		33x18x29cm		36x27x31cm		57x33x33cm	
ULTRASONIC POWER	100W		180W		200W		500W	
HEATING POWER	-	88W	-	170W	-	270W	-	550W

### Ordering Information

PRODUCT CODE	DESCRIPTION	PACK SIZE
<a href="#">188203-22</a>	3L Ultrasonic Bath	Each
<a href="#">188204-22</a>	3L Ultrasonic Bath with Heating	Each
<a href="#">188206-22</a>	6L Ultrasonic Bath	Each
<a href="#">188207-22</a>	6L Ultrasonic Bath with Heating	Each
<a href="#">188210-22</a>	10L Ultrasonic Bath	Each
<a href="#">188211-22</a>	10L Ultrasonic Bath with Heating	Each
<a href="#">188220-22</a>	20L Ultrasonic Bath	Each
<a href="#">188221-22</a>	20L Ultrasonic Bath with Heating	Each



# Accessories

## WIRE BASKET

Ideal for smaller parts.

### Ordering Information

PRODUCT CODE	DESCRIPTION	PACK SIZE
<a href="#">188203-30</a>	Wire Basket for 3L	Each
<a href="#">188206-30</a>	Wire Basket for 6L	Each
<a href="#">188210-30</a>	Wire Basket for 10L	Each
<a href="#">188220-30</a>	Wire Basket for 20L	Each




## WIRE BRACKET

Ideal for larger parts.

### Ordering Information

PRODUCT CODE	DESCRIPTION	PACK SIZE
<a href="#">188203-31</a>	Wire Bracket for 3L	Each
<a href="#">188206-31</a>	Wire Bracket for 6L	Each
<a href="#">188210-31</a>	Wire Bracket for 10L	Each
<a href="#">188220-31</a>	Wire Bracket for 20L	Each

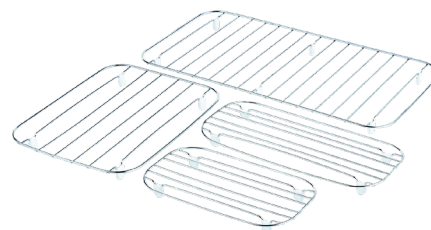

## BEAKER COVER

Used to position beaker for indirect cleaning.

### Ordering Information

PRODUCT CODE	DESCRIPTION	PACK SIZE
<a href="#">188200-31</a>	Beaker Cover for 3L and 6L	Each
<a href="#">188200-32</a>	Beaker Cover for 10L	Each
<a href="#">188200-33</a>	Beaker Cover for 20L	Each





Printed technical specifications are subject to change without notification.  
LabCo® is a registered trademark.

See the full product range online.

1800 052 226  
sales@labcoscientific.com.au  
[labcoscientific.com.au](http://labcoscientific.com.au)